

Surgical protocol for copaSKY Ø 3.0



Implant length

8, 10, 12, 14 mm

Crestal drill

300 rpm

REF COPACD30

3.0

Hard bone D1

300 rpm

3.0

Medium hard /
soft bone D2-D4

300 rpm

3.0

Pilot drill

Twistdrill

REF SKY-DP06

REF SKY-DP08

REF SKYDT23K

REF SKYDT23L

800 - 1.000 rpm

Surgical protocol for copaSKY

Implant length

8, 10
12, 14 mm

5.2, 8, 10
12, 14 mm

5.2, 8, 10
12, 14 mm

5.2
8 mm

5.2
8 mm

Crestal drill

300 rpm

REF SKYCD35n

3.5 N

REF SKYCD40

4.0

REF SKYCD45

4.5

REF SKYCD55

5.0

REF COPACD60

6.0

Hard bone D1

300 rpm

3.5 N

4.0

4.5

5.0

6.0

Medium hard /
soft bone D2-D4

300 rpm

3.5 N

4.0

4.5

5.0

6.0

Pilot drill

Twistdrill

REF SKY-DP06

REF SKY-DP08

REF SKYDT23K

REF SKYDT23L

800 - 1.000 rpm

Attention:
The laser marking of the
pilot drill indicates the
drill depth for the short
copaSKY implants

Caution!
The implant axis may
deviate slightly from the
cavity axis.

bredent
medical

bredent
group

bredent
medical

bredent
group



Surgical protocol for

narrow
SKY

blue
SKY

classic
SKY

white
SKY

- optimized, bone quality-oriented set of instruments and surgical protocol for unsurpassed primary stability
- Drills with detachable drill stops
- Reduction of the number of drills for increased control and dependability during surgery

Implant length

8, 10, 12, 14, 16 mm

8, 10, 12, 14, 16 mm

8, 10, 12, 14, mm

8, 10, 12, mm



Crestal drill

300 rpm

REF SKYCD35n

REF SKYCD40

REF SKYCD45

REF SKYCD55

Hard bone D1

300 rpm

REF SKYD1235

REF SKYD1240

REF SKYD1245

REF SKYD1255

Medium hard / soft bone D2-D4

300 rpm

REF SKYD3435

REF SKYD3440

REF SKYD3445

REF SKYD3455

Pilot drill

Twistdrill

REF SKY-DP06

REF SKY-DP08


REF SKYD123K

REF SKYD123L


800 - 1.000 rpm

With the 3.5 SKY classic implants, the crestal drill is only sunk up to the laser marking.

In the following implants:
• narrowSKY
• blueSKY 4.0 to 5.5
• SKY classic 4.0 to 4.5
the crestal drill is completely inserted.

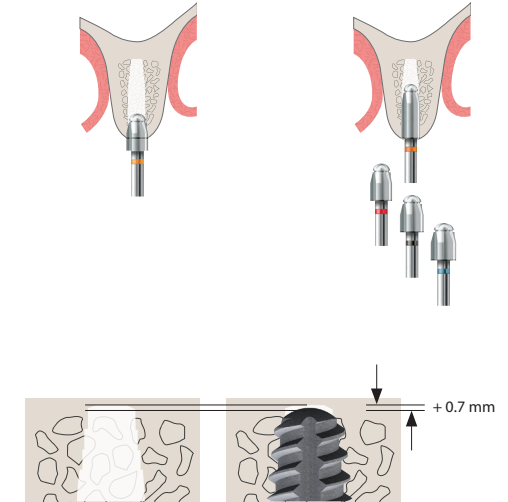

Detachable drill stops

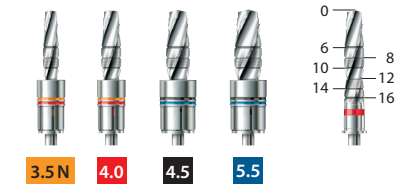

bredent
medical

bredent
dental